



01 Digital Abutment


Evolved from stock abutment, simple yet, compensate the weak points of Customized Abutment It's Digital Abutment

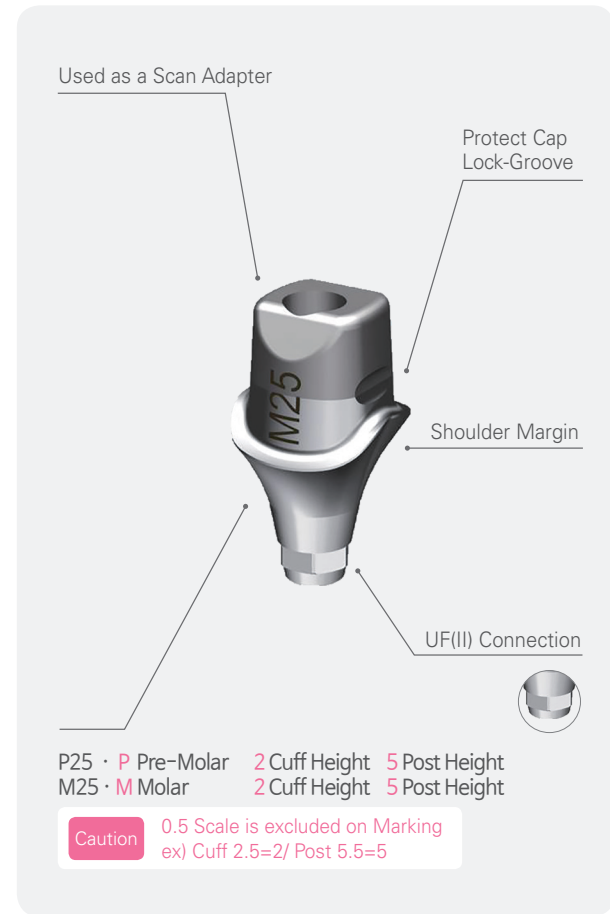

02 Connecting Hex becomes simple and easy

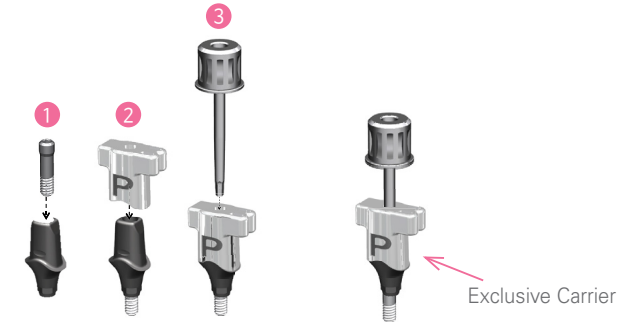

Easier delivery to the intra oral area (Only require the driver)

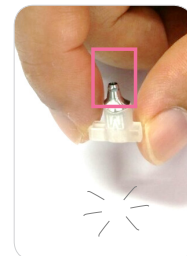


The Carrier prevents screw from dropping into the patient's throat.



It is designed for easier handling and precise connecting.

03 Digital Prosthesis Process

Full Digital Abutment utilizes the detailed 3D CAD system and it changes whole process of making prosthesis in analogue method.

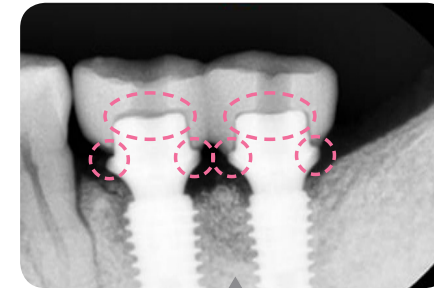

Analogue method

If the margin line is subgingival, it is hard to take accurate impression and it leads to inaccuracy of prosthesis.

VS



Digital Process

Even if the margin line is subgingival, digital process uses a library system to make an ideal result. It doesn't even require the gingiva retraction cord.

04 Healing a gingiva and making prosthesis at the same time



» Exclusive Healing Cap for DA makes it easier to recover the gingiva and to scan intraoral site.



» Exclusive Healing Cap for DA helps to keep the margin line during the healing stage of gingiva and it needs no packing cord.




Final Prosthesis

TIP
DA contains functions of healing abutment, scan body and the abutment for final prosthesis. It means DA is economical and also decrease the chair-time.

Can't be better, It is Digital Abutment

The ideal size of DA for the each patients will be sent to you prior to the surgery. This is also applied to DIONavi. surgical guide.



Making a prosthesis using Digital Abutment & DIONavi.

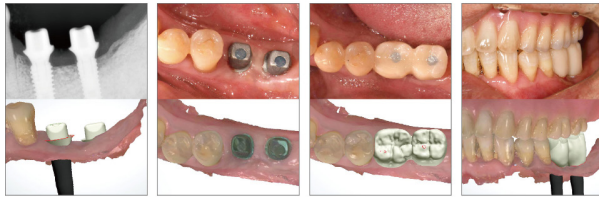

- It is possible to identify the exact height of gingiva during the surgical planning process and it enables the right selection of DA prior to the surgery.
- Surgical guide, DA, Abutment Jig for the immediate restoration case can be made in advance of the surgery.
- Even if it needs to be replaced to new cuff in another size because of changed gingival area, the patient still doesn't need to visit the clinic to take new impression but you can choose proper sized DA by using a library system.


Pre-Molar

No.	Cuff	Post
1	1.5	4.0
2	2.5	
3	3.5	
4	4.5	
5	1.5	5.5
6	2.5	
7	3.5	
8	4.5	

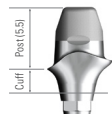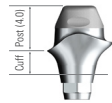

Molar

No.	Cuff	Post
1	1.5	4.0
2	2.5	
3	3.5	
4	4.5	
5	1.5	5.5
6	2.5	
7	3.5	
8	4.5	

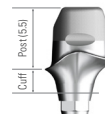

DIO Implant Online Community

www.facebook.com/diohq

www.youtube.com/DIOHQ

DIGITAL ABUTMENT

What makes Digital Abutment so special?

Advantages of Stock Abutment

- Shorter making process
- Lower cost
- Well-fitted prosthesis by using library ...system

Advantages of Customized Abutment

- Better protect the prosthesis from rotation than rounded stock abutment.
- Customized for the patient's gingiva
- Lower chance of getting omitted prosthesis due to the wider post area

Easy to use like stock abutment,
Customized for the individual like
custom abutment

All in One!
Digital Abutment



DIO IMPLANT www.dioimplant.com
DIO Implant, 66 Centum Seo-ro, Haeundae-gu, Busan
T +82 51 745 7777 F +82 51 745 7778