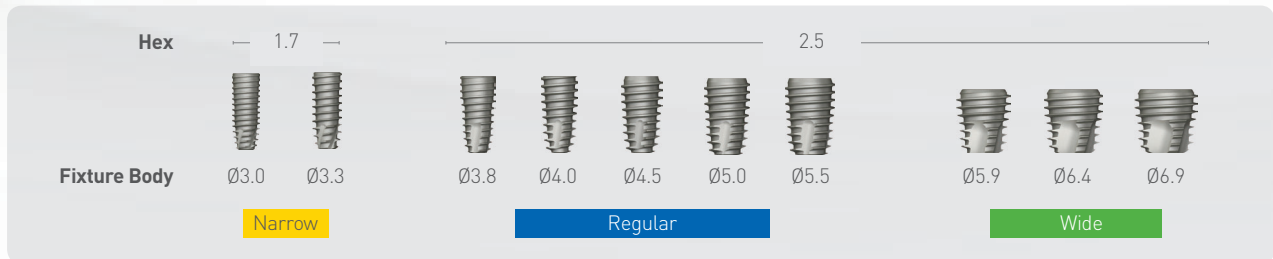


UFII *IMPLANT*

Internal Submerged System

DIO IMPLANT

UFII Int. Submerged System

UFII Implant Advantage**Open Thread**

Fixture can be placed with minimal resistance, allowing depth adjustment

Platform Switching

Platform Switching structure allows stable and safe attachment

Upper Thread**Taper Thread, Initial Fixation**

Taper threaded straight body allows effective stress-distribution. Friction can be reduced in the cortical bone with the shallow root depth design going upwards

Morse Tapered Contact Connection

Provides an excellent level of connectivity and loading distribution of fixture and abutment

Straight Body**Minimizes Bone Damage**

Fixture length changes based on body length hence, long fixtures allow stable placement in the same level as straight types

Apical Design**Tapered Apical Design**


Initial insertion can be executed easily even for weak bone conditions, with the tapered design

Cutting edge

Sharp cutting edge increases self-tapping efficiency

Apex

Rounded apex design minimizes damage to the inferior alveolar nerve and sinus floor

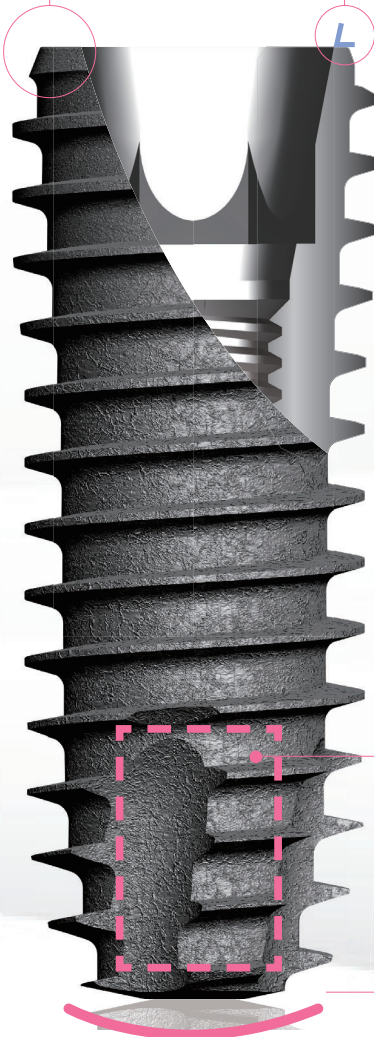

UFII Int. Submerged System

UFII Implant System

UFII Narrow

D Ø3.0 | Ø3.3

L 8.5 | 10 | 11.5 | 13 | 15

HSA

UV Active

UFII Regular

D Ø3.8 | Ø4.0 | Ø4.5 | Ø5.0 | Ø5.5

L 7 | 8.5 | 10 | 11.5 | 13 | 15 | 16 | 18

HSA

UV Active


UFII Wide

D Ø5.9 | Ø6.4 | Ø6.9

L 7 | 8.5 | 10 | 11.5 | 13 | 15


HSA

UV Active

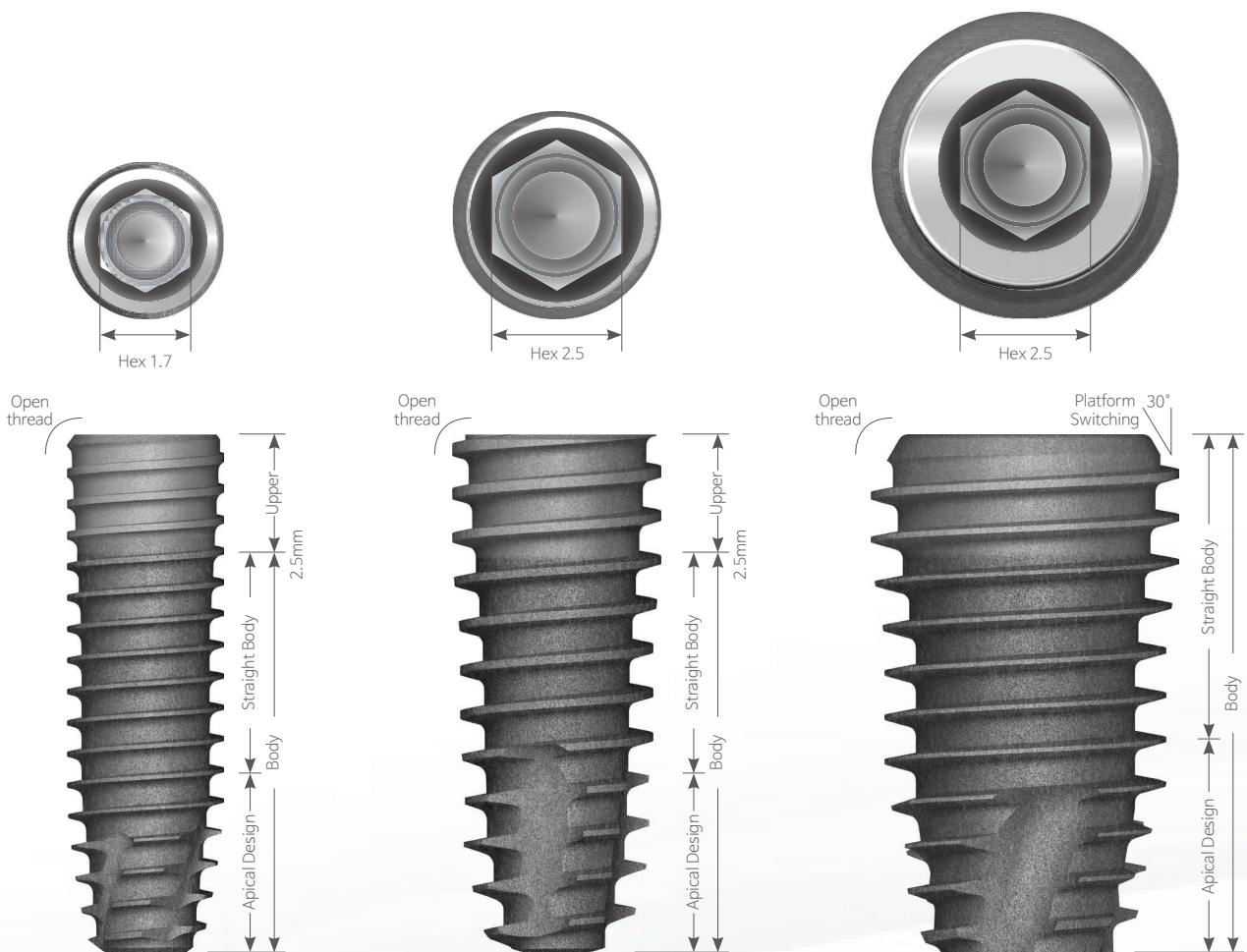

Hybrid Type SLA Implant

UFII Implant Surface **HSA**

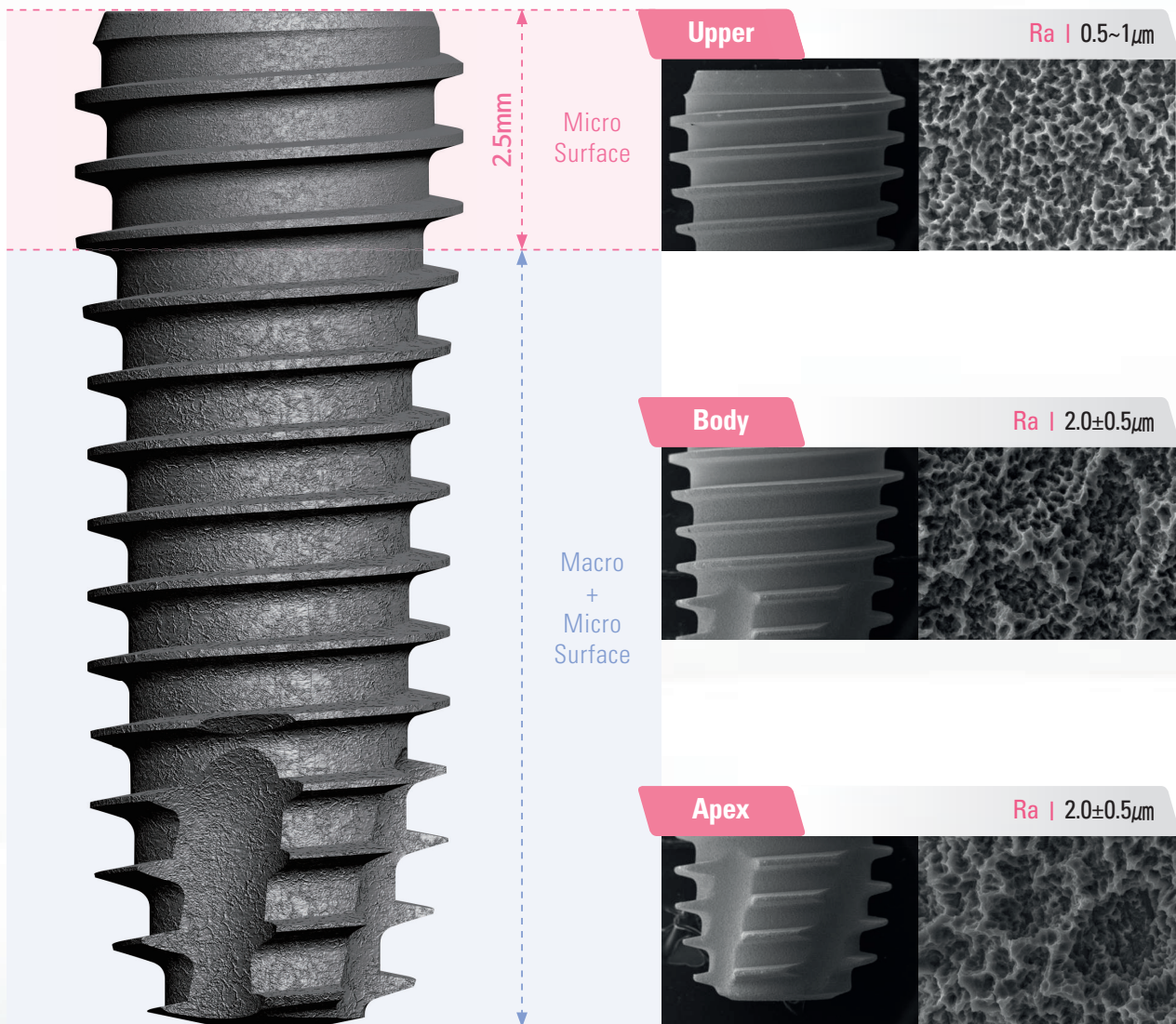

IMPLANT UFII SYSTEM


Narrow

Diameter : Ø3.0 / Ø3.3

Length : 8.5 / 10 / 11.5 / 13 / 15


Ø3.0 Fixture : Only for anterior

Implant Surface : HSA / UV Active

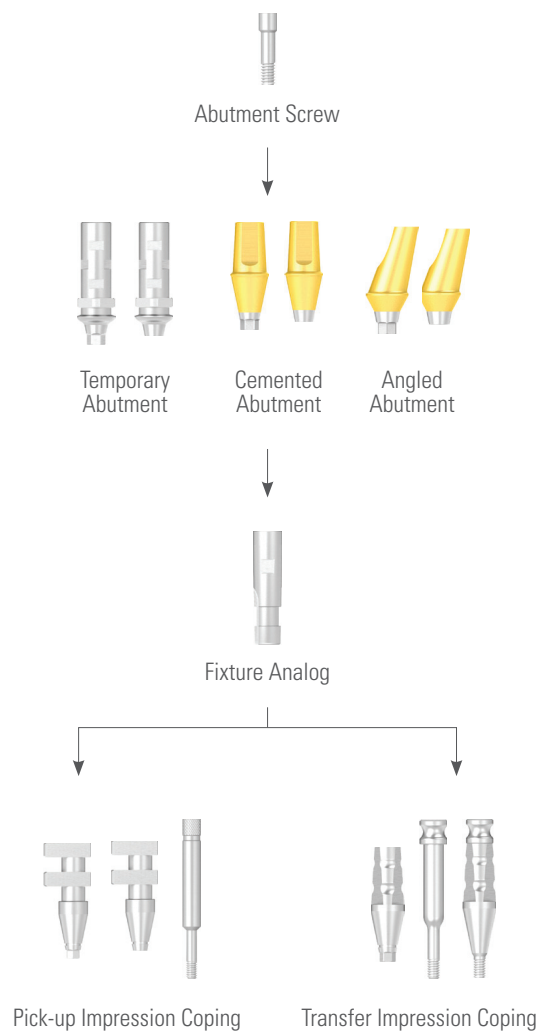

Cement-Retained Restorations

Temporary / Cemented / Angled Abutment

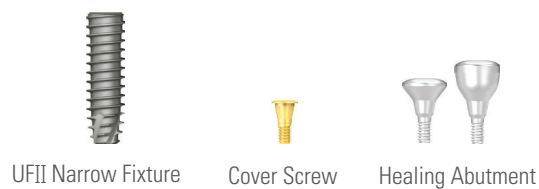



Single Bridge
Plastic Coping



Solid Abutment Analog



Impression Cap



Protect Cap

Solid Abutment

Cement-Retained Restorations

Solid Abutment



Stainless Steel Housing



OT Cap



Ball Abutment Analog



Ball Abutment

Overdenture-Retained Restorations

Ball Abutment



UFII Narrow Fixture



Cover Screw



Healing Abutment

Narrow Fixture

- HSA / Packing Unit :
 Fixture | Order code : Product Code + SF (ex : UF(II)N 3008SF)
 Fixture & Coverscrew | Order code : Product Code + S (ex : UF(II)N 3008S)
- UV Active / Packing Unit : Fixture

HSA

Hex	Length	Ø3.0	Ø3.3
1.7	8.5	UF(II)N 3008S	UF(II)N 3308S
	10	UF(II)N 3010S	UF(II)N 3310S
	11.5	UF(II)N 3011S	UF(II)N 3311S
	13	UF(II)N 3013S	UF(II)N 3313S
	15	UF(II)N 3015S	UF(II)N 3315S
Apex		Ø1.6	Ø1.9

UV Active

Hex	Length	Ø3.0	Ø3.3
1.7	8.5	UF(II)N 3008UV3	UF(II)N 3308UV3
	10	UF(II)N 3010UV3	UF(II)N 3310UV3
	11.5	UF(II)N 3011UV3	UF(II)N 3311UV3
	13	UF(II)N 3013UV3	UF(II)N 3313UV3
	15	-	-
Apex		Ø1.6	Ø1.9

Healing Abutment

Height	Cuff	Ø4.1	Ø4.2	Ø4.6	Ø4.7
2	0	UNSHA 4020		UNSHA 4520	
4	2		UNSHA 4024		UNSHA 4524
5.5	3		UNSHA 4035		UNSHA 4535
7	4		UNSHA 4047		UNSHA 4547
8.5	5		UNSHA 4058		UNSHA 4558
10	6		UNSHA 4060		UNSHA 4560

Temporary Abutment

Cuff	Ø4.0	
	Hex (Guide Pin)	Non-Hex
1	UNSTA 4010H	UNSTA 4010N
3	UNSTA 4030H	UNSTA 4030N

Cemented Abutment

Cuff	Ø4.0			
	Hex		Non-Hex	
	L 5.5	L 7	L 5.5	L 7
1	UNSCA 40105H	UNSCA 40107H	UNSCA 40105N	UNSCA 40107N
2	UNSCA 40205H	UNSCA 40207H	UNSCA 40205N	UNSCA 40207N
3	UNSCA 40305H	UNSCA 40307H	UNSCA 40305N	UNSCA 40307N
4	UNSCA 40405H	UNSCA 40407H	UNSCA 40405N	UNSCA 40407N
5	UNSCA 40505H	UNSCA 40507H	UNSCA 40505N	UNSCA 40507N
6	UNSCA 40605H	UNSCA 40607H	UNSCA 40605N	UNSCA 40607N
7	UNSCA 40705H	UNSCA 40707H	UNSCA 40705N	UNSCA 40707N

Angled Abutment

Cuff	Ø4.0		
	Hex A	Hex B	Non-Hex
2	UNSAA 402015AH	UNSAA 402015BH	UNSAA 402015N
4	UNSAA 404015AH	UNSAA 404015BH	UNSAA 404015N

Solid Abutment

Cuff	Ø4.0		
	L 4.0	L 5.5	L 7
1	UNSSA 40104	UNSSA 40105	UNSSA 40107
2	UNSSA 40204	UNSSA 40205	UNSSA 40207
3	UNSSA 40304	UNSSA 40305	UNSSA 40307
4	UNSSA 40404	UNSSA 40405	UNSSA 40407
5	UNSSA 40504	UNSSA 40505	UNSSA 40507
6	UNSSA 40604	UNSSA 40605	UNSSA 40607
7	UNSSA 40704	UNSSA 40705	UNSSA 40707

Ball Abutment

Length	Ø4.0
1	UNSBA 3510
2	UNSBA 3520
3	UNSBA 3530
4	UNSBA 3540
5	UNSBA 3550
6	UNSBA 3560
7	UNSBA 3570

IMPLANT UFII SYSTEM

Regular

Diameter : Ø3.8 / Ø4.0 / Ø4.5 / Ø5.0 / Ø5.5

Length : 7 / 8.5 / 10 / 11.5 / 13

Implant Surface : HSA / UV Active


IMPLANT UFII SYSTEM

Wide

Diameter : Ø5.9 / Ø6.4 / Ø6.9

Length : 7 / 8.5 / 10 / 11.5


Implant Surface : HSA / UV Active

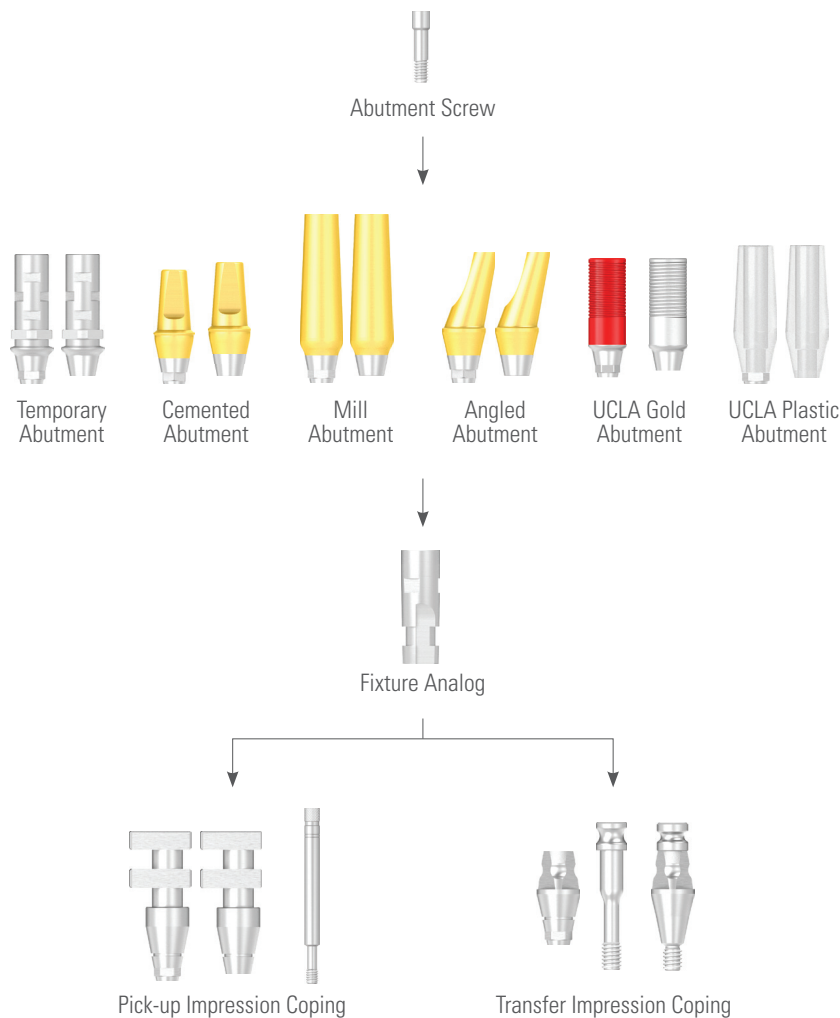
**Cement-Retained Restorations**

Temporary / Cemented / Mill / Angled Abutment

Screw-Retained Restorations


UCLA Gold / UCLA Plastic Abutment

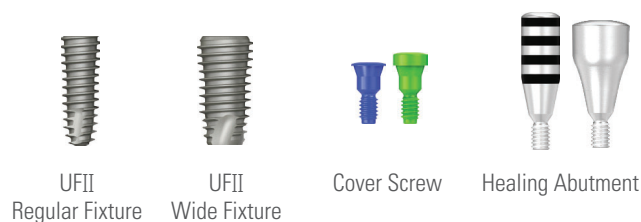

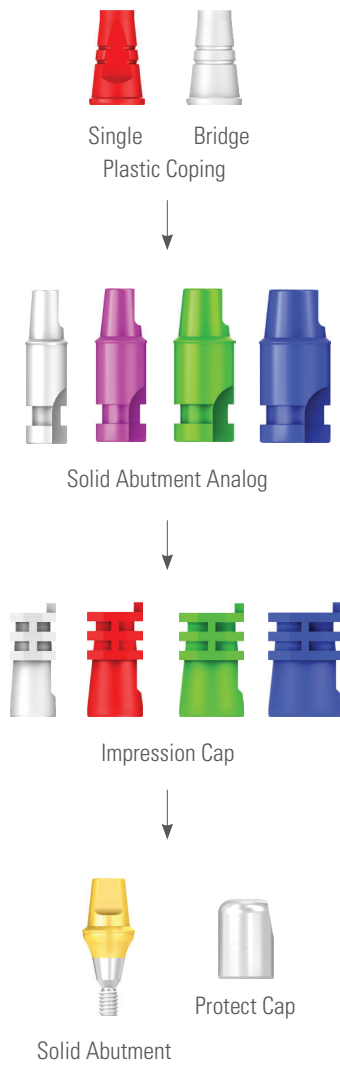


Cement-Retained Restorations



Solid Abutment

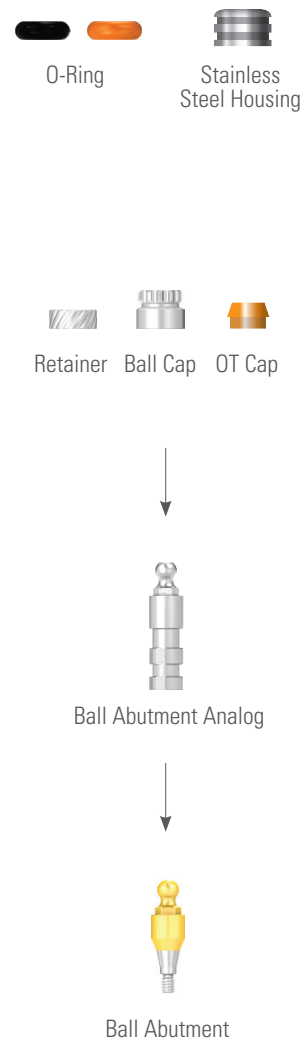

Overdenture-Retained Restorations

Ball Abutment

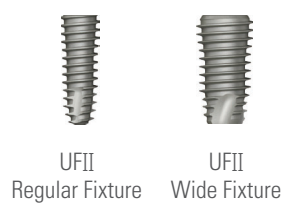

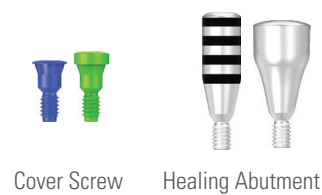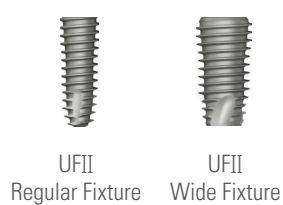

Regular / Wide Fixture

- HSA / Packing Unit :
Fixture | Order code : Product Code + SF (ex : UF(II) 3808SF)
Fixture & Coverscrew | Order code : Product Code + S (ex : UF(II) 3808S)
- UV Active / Packing Unit : Fixture

HSA Regular

Hex	Length	Ø3.8	Ø4.0	Ø4.5	Ø5.0	Ø5.5
2.5	7	-	-	UF(II) 4507S	UF(II) 5007S	UF(II) 5507S
	8.5	UF(II) 3808S	UF(II) 4008S	UF(II) 4508S	UF(II) 5008S	UF(II) 5508S
	10	UF(II) 3810S	UF(II) 4010S	UF(II) 4510S	UF(II) 5010S	UF(II) 5510S
	11.5	UF(II) 3811S	UF(II) 4011S	UF(II) 4511S	UF(II) 5011S	UF(II) 5511S
	13	UF(II) 3813S	UF(II) 4013S	UF(II) 4513S	UF(II) 5013S	UF(II) 5513S
Apex		Ø2.1	Ø2.2	Ø2.6	Ø3.1	Ø3.7

UV Active Regular

Hex	Length	Ø3.8	Ø4.0	Ø4.5	Ø5.0	Ø5.5
2.5	7	-	-	UF(II) 4507UV3	UF(II) 5007UV3	UF(II) 5507UV3
	8.5	UF(II) 3808UV3	UF(II) 4008UV3	UF(II) 4508UV3	UF(II) 5008UV3	UF(II) 5508UV3
	10	UF(II) 3810UV3	UF(II) 4010UV3	UF(II) 4510UV3	UF(II) 5010UV3	UF(II) 5510UV3
	11.5	UF(II) 3811UV3	UF(II) 4011UV3	UF(II) 4511UV3	UF(II) 5011UV3	UF(II) 5511UV3
	13	UF(II) 3813UV3	UF(II) 4013UV3	UF(II) 4513UV3	UF(II) 5013UV3	UF(II) 5513UV3
Apex		Ø2.1	Ø2.2	Ø2.6	Ø3.1	Ø3.7

HSA Wide

Hex	Length	Ø5.9	Ø6.4	Ø6.9
2.5	7	UF(II) 6007S	UF(II) 6507S	UF(II) 7007S
	8.5	UF(II) 6008S	UF(II) 6508S	UF(II) 7008S
	10	UF(II) 6010S	UF(II) 6510S	UF(II) 7010S
	11.5	UF(II) 6011S	UF(II) 6511S	UF(II) 7011S
Apex		Ø5.4	Ø5.9	Ø6.4

UV Active Wide

Hex	Length	Ø5.9	Ø6.4	Ø6.9
2.5	7	UF(II) 6007UV3	UF(II) 6507UV3	-
	8.5	UF(II) 6008UV3	UF(II) 6508UV3	-
	10	UF(II) 6010UV3	UF(II) 6510UV3	-
	11.5	UF(II) 6011UV3	UF(II) 6511UV3	-
Apex		Ø5.4	Ø5.9	Ø6.4

Healing Abutment

Cuff	Ø4.0
2	SSHA 4020
4	SSHA 4040
6	SSHA 4060
8	SSHA 4080

Height	Cuff	D1 - Ø4.6 D2 - Ø4.5	D1 - Ø5.6 D2 - Ø5.5	D1 - Ø6.6 D2 - Ø6.5	D1 - Ø7.6 D2 - Ø7.5
2	2	-	-	-	-
4	2	SSHA 4524	SSHA 5524	SSHA 6524	SSHA 7524
5.5	3	SSHA 4535	SSHA 5535	SSHA 6535	SSHA 7535
7	4	SSHA 4547	SSHA 5547	SSHA 6547	SSHA 7547

Temporary Abutment

Cuff	Ø4.5
Hex	Non-Hex
1	SSTA 4510H[H] SSTA 4510N[H]
3	SSTA 4530H[H] SSTA 4530N[H]

Mill Abutment

Cuff	Diameter	Ø4.5	
		Hex	Non-Hex
1.5	Ø4.0	SSMA 4015H[H]	SSMA 4015N[H]
2	Ø4.5	SSMA 4520H[H]	SSMA 4520N[H]
2.5	Ø5.5	SSMA 5525H[H]	SSMA 5525N[H]
3	Ø6.5	SSMA 6530H[H]	SSMA 6530N[H]
3	Ø7.5	SSMA 7530H[H]	SSMA 7530N[H]

UCLA Abutment

Gold Type

Cuff	Ø4.5
Hex	Non Hex
1	SSGA 4510H[H] SSGA 4510N[H]
3	SSGA 4530H[H] SSGA 4530N[H]

UCLA Plastic Abutment

Cuff	Ø4.5		Ø5.5	
	Hex	Non Hex	Hex	Non Hex
1	SSPA 4512H	SSPA 4512N	SSPA 5512H	SSPA 5512N

Octa Abutment

Cuff	Ø4.8	Ø6.0	Ø6.5
1.5	SSOA 4815	SSOA 6015	SSOA 6515
2.5	SSOA 4825	SSOA 6025	SSOA 6525
3.5	SSOA 4835	SSOA 6035	SSOA 6535
4.5	SSOA 4845	SSOA 6045	SSOA 6545
5.5	SSOA 4855	SSOA 6055	SSOA 6555

Ball Abutment

Length	Ø3.5
1	SSBA 3510
2	SSBA 3520
3	SSBA 3530
4	SSBA 3540
5	SSBA 3550
6	SSBA 3560

Multiunit Abutment

Multiunit Straight Abutment

Length	Ø4.8
1.5	MSA 4801
2.5	MSA 4802
3.5	MSA 4803
4.5	MSA 4804

Multiunit Angle Abutment

Cuff	Angle 20°	Angle 30°
2.5	MAA 482520H	-
3	MAA 483520H	MAA 483530H
4	MAA 484520H	MAA 484530H
5	MAA 485520H	MAA 485530H

Angled Abutment

		Ø4.5				
Angled	Type	1.5	2	3	4	5
15°	Hex A	SSAA 451515AH[H]	SSAA 452015AH[H]	SSAA 453015AH[H]	SSAA 454015AH[H]	SSAA 455015AH[H]
	Hex B	SSAA 451515BH[H]	SSAA 452015BH[H]	SSAA 453015BH[H]	SSAA 454015BH[H]	SSAA 455015BH[H]
	Non-Hex	SSAA 451515N[H]	SSAA 452015N[H]	SSAA 453015N[H]	SSAA 454015N[H]	SSAA 455015N[H]
25°	Hex A	SSAA 451525AH[H]	SSAA 452025AH[H]	SSAA 453025AH[H]	SSAA 454025AH[H]	SSAA 455025AH[H]
	Hex B	SSAA 451525BH[H]	SSAA 452025BH[H]	SSAA 453025BH[H]	SSAA 454025BH[H]	SSAA 455025BH[H]
	Non-Hex	SSAA 451525N[H]	SSAA 452025N[H]	SSAA 453025N[H]	SSAA 454025N[H]	SSAA 455025N[H]

		Ø5.5				
Angled	Type	1.5	2	3	4	5
15°	Hex A	SSAA 551515AH[H]	SSAA 552015AH[H]	SSAA 553015AH[H]	SSAA 554015AH[H]	SSAA 555015AH[H]
	Hex B	SSAA 551515BH[H]	SSAA 552015BH[H]	SSAA 553015BH[H]	SSAA 554015BH[H]	SSAA 555015BH[H]
	Non-Hex	SSAA 551515N[H]	SSAA 552015N[H]	SSAA 553015N[H]	SSAA 554015N[H]	SSAA 555015N[H]
25°	Hex A	SSAA 551525AH[H]	SSAA 552025AH[H]	SSAA 553025AH[H]	SSAA 554025AH[H]	SSAA 555025AH[H]
	Hex B	SSAA 551525BH[H]	SSAA 552025BH[H]	SSAA 553025BH[H]	SSAA 554025BH[H]	SSAA 555025BH[H]
	Non-Hex	SSAA 551525N[H]	SSAA 552025N[H]	SSAA 553025N[H]	SSAA 554025N[H]	SSAA 555025N[H]

Solid Abutment

Length	Cuff	Ø4.5	Ø5.5	Ø6.5	Ø7.5
4	1	SSSA 45104	SSSA 55104	SSSA 65104	SSSA 75104
	2	SSSA 45204	SSSA 55204	SSSA 65204	SSSA 75204
	3	SSSA 45304	SSSA 55304	SSSA 65304	SSSA 75304
	4	SSSA 45404	SSSA 55404	SSSA 65404	SSSA 75404
	5	SSSA 45504	SSSA 55504	SSSA 65504	SSSA 75504
5.5	1	SSSA 45105	SSSA 55105	SSSA 65105	SSSA 75105
	2	SSSA 45205	SSSA 55205	SSSA 65205	SSSA 75205
	3	SSSA 45305	SSSA 55305	SSSA 65305	SSSA 75305
	4	SSSA 45405	SSSA 55405	SSSA 65405	SSSA 75405
	5	SSSA 45505	SSSA 55505	SSSA 65505	SSSA 75505
7	1	SSSA 45107	SSSA 55107	SSSA 65107	SSSA 75107
	2	SSSA 45207	SSSA 55207	SSSA 65207	SSSA 75207
	3	SSSA 45307	SSSA 55307	SSSA 65307	SSSA 75307
	4	SSSA 45407	SSSA 55407	SSSA 65407	SSSA 75407
	5	SSSA 45507	SSSA 55507	SSSA 65507	SSSA 75507

Cemented Abutment

Hex Type

Length	Cuff	Ø4.5	Ø5.5	Ø6.5	Ø7.5
4	1	SSCA 45104H(II)[H]	SSCA 55104H(III)[H]	SSCA 65104H(III)[H]	SSCA 75104H(III)[H]
	2	SSCA 45204H(II)[H]	SSCA 55204H(III)[H]	SSCA 65204H(III)[H]	SSCA 75204H(III)[H]
	3	SSCA 45304H(II)[H]	SSCA 55304H(III)[H]	SSCA 65304H(III)[H]	SSCA 75304H(III)[H]
	4	SSCA 45404H(II)[H]	SSCA 55404H(III)[H]	SSCA 65404H(III)[H]	SSCA 75404H(III)[H]
	5	SSCA 45504H(II)[H]	SSCA 55504H(III)[H]	SSCA 65504H(III)[H]	SSCA 75504H(III)[H]
5.5	1	SSCA 45105H(II)[H]	SSCA 55105H(III)[H]	SSCA 65105H(III)[H]	SSCA 75105H(III)[H]
	2	SSCA 45205H(II)[H]	SSCA 55205H(III)[H]	SSCA 65205H(III)[H]	SSCA 75205H(III)[H]
	3	SSCA 45305H(II)[H]	SSCA 55305H(III)[H]	SSCA 65305H(III)[H]	SSCA 75305H(III)[H]
	4	SSCA 45405H(II)[H]	SSCA 55405H(III)[H]	SSCA 65405H(III)[H]	SSCA 75405H(III)[H]
	5	SSCA 45505H(II)[H]	SSCA 55505H(III)[H]	SSCA 65505H(III)[H]	SSCA 75505H(III)[H]
7	1	SSCA 45107H(II)[H]	SSCA 55107H(III)[H]	SSCA 65107H(III)[H]	SSCA 75107H(III)[H]
	2	SSCA 45207H(II)[H]	SSCA 55207H(III)[H]	SSCA 65207H(III)[H]	SSCA 75207H(III)[H]
	3	SSCA 45307H(II)[H]	SSCA 55307H(III)[H]	SSCA 65307H(III)[H]	SSCA 75307H(III)[H]
	4	SSCA 45407H(II)[H]	SSCA 55407H(III)[H]	SSCA 65407H(III)[H]	SSCA 75407H(III)[H]
	5	SSCA 45507H(II)[H]	SSCA 55507H(III)[H]	SSCA 65507H(III)[H]	SSCA 75507H(III)[H]
Screw		SSC 2008SH	SSC 2008H	SSC 2008H	SSC 2008H

Non-Hex Type

Length	Cuff	Ø4.5	Ø5.5	Ø6.5	Ø7.5
4	1	SSCA 45104N(II)[H]	SSCA 55104N(III)[H]	SSCA 65104N(III)[H]	SSCA 75104N(III)[H]
	2	SSCA 45204N(II)[H]	SSCA 55204N(III)[H]	SSCA 65204N(III)[H]	SSCA 75204N(III)[H]
	3	SSCA 45304N(II)[H]	SSCA 55304N(III)[H]	SSCA 65304N(III)[H]	SSCA 75304N(III)[H]
	4	SSCA 45404N(II)[H]	SSCA 55404N(III)[H]	SSCA 65404N(III)[H]	SSCA 75404N(III)[H]
	5	SSCA 45504N(II)[H]	SSCA 55504N(III)[H]	SSCA 65504N(III)[H]	SSCA 75504N(III)[H]
5.5	1	SSCA 45105N(II)[H]	SSCA 55105N(III)[H]	SSCA 65105N(III)[H]	SSCA 75105N(III)[H]
	2	SSCA 45205N(II)[H]	SSCA 55205N(III)[H]	SSCA 65205N(III)[H]	SSCA 75205N(III)[H]
	3	SSCA 45305N(II)[H]	SSCA 55305N(III)[H]	SSCA 65305N(III)[H]	SSCA 75305N(III)[H]
	4	SSCA 45405N(II)[H]	SSCA 55405N(III)[H]	SSCA 65405N(III)[H]	SSCA 75405N(III)[H]
	5	SSCA 45505N(II)[H]	SSCA 55505N(III)[H]	SSCA 65505N(III)[H]	SSCA 75505N(III)[H]

Specifications & Implant Size

Recommendations

Tooth Number	#9	#10	#11	#12	#13	#14	#15
Fixture Diameter	Ø3.8	Ø3.8	Ø4.5	Ø4.5	Ø4.5	Ø5.5	Ø5.5
Abut. Diameter	Ø4.5	Ø4.5	Ø4.5	Ø5.5	Ø5.5	Ø6.5	Ø6.5



Tooth Number	#24	#23	#22	#21	#20	#19	#18
Fixture Diameter	Ø3.0	Ø3.3	Ø4.0	Ø4.5	Ø5.0	Ø5.5	Ø5.5
Abut. Diameter	Ø4.0	Ø4.0	Ø4.5	Ø5.5	Ø5.5	Ø6.5	Ø6.5

Placement after molar extraction

Fixture Diameter	Ø5.9	Ø6.4	Ø6.9
------------------	------	------	------



Insufficient vertical bone quantity


Fixture Diameter	Ø4.5	Ø5.0	Ø5.5	Ø5.9	Ø6.4	Ø6.9
Fixture Length	7mm					



UFII Int. Submerged System Surgical Kit

Guided Surgical Kits

DIONavi. Master Kit

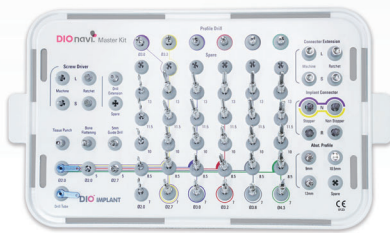

DIONavi. Narrow Kit




DIONavi. Wide Kit



DIONavi. Flapless Crestal Sinus Kit




DIONavi. Special Kit

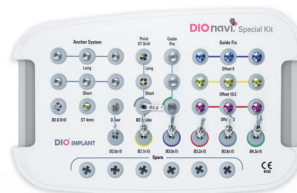


Non-Guided Surgical Kits


UF II Master kit

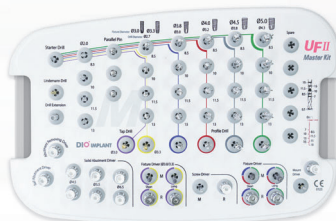


UF II Surgical Kit

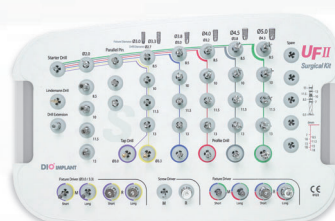


UF II Narrow Kit

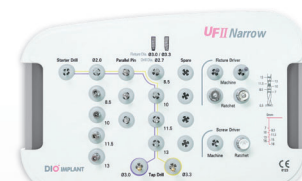


UF II Wide Kit

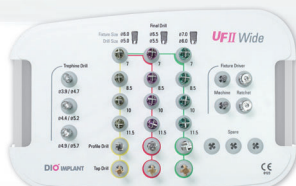

Dr. SOS+ Kit

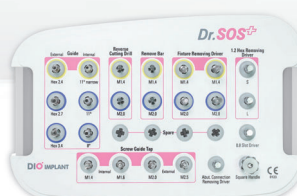

UF II Prosthetic Kit

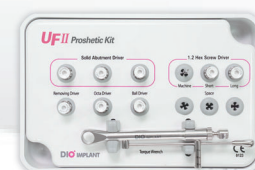



www.dioimplant.com

DIO USA

3470 Wilshire blvd.,#620,

Los Angeles

CA 90010 USA

T 213-365-2875

E usahq@dioimplantusa.com